


Now Available From Hummingbird Med

Disposable Closed Blood Draw and Flush Set for Central Lines

The **SAMMI SET**.... **S**ingle **A**ccess to any CL catheter for blood **M**anagement and fluid **M**anagement providing **I**mproved line care and **S**tandardization for **E**very blood **T**ransfer procedure.

BTF-HM-3 SAMMI Set (Blood Transfer and Flush Set)

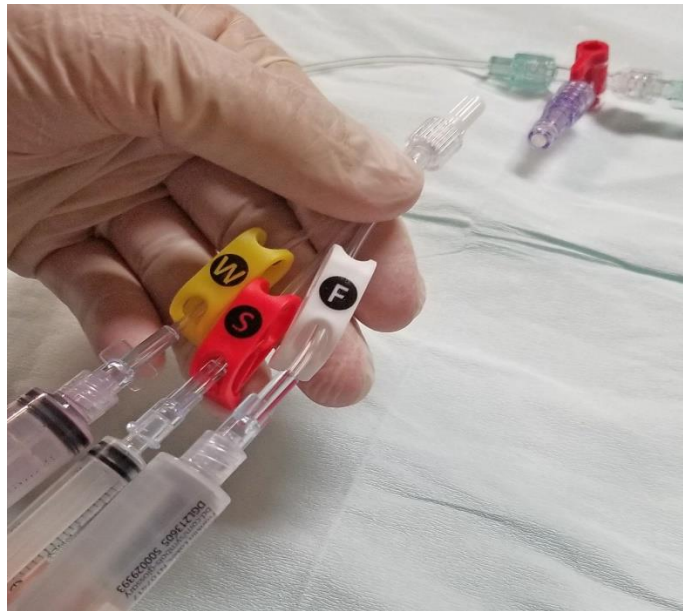

- The SAMMI set is attached to a single (closed) port access point in the vascular access line or catheter anytime or location that a blood draw is needed from a Central Line catheter. (ICU Med Surg, ER, Outpatient areas)
- Blood Draw is completed in a closed procedure at the access point with only one access to the line or catheter for Pre-Flush, Waste, Sample and Post Flush. (pre-flush is optional)
- Improves line maintenance and standardizes blood transfer and flush for every blood draw.
- Never open to atmosphere during use, providing fully closed system in use. One time use.
- Provides for blood conservation in use as residual blood discard is minimal. (less than 0.3mL)
- Return of Waste helps minimize anemia risk in critical patients. (waste return is optional)
- Can be used in any circumstance where a Closed Blood Draw is needed from a central line and access is available. (ie. KVO, Hep-Locked Catheters) Single Use Disposable.

Order Number: BTF-HM-3 50 per box